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Back and neck warmer


with massage function

Feel the pleasant warmth, for example


... at work

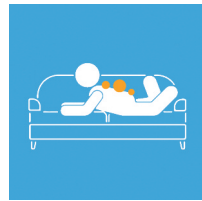... on the couch

... in the kitchen

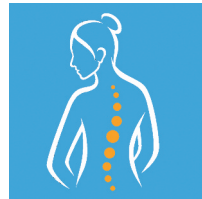with massage function

Improve your blood circulation through the gentle warmth and the relaxing massage function!

Well-functioning blood circulation is essential for a fit and healthy body. The body must constantly supply its organs and tissue with oxygen and nutrients and transport away metabolic waste again. Good blood circulation is important for the full functioning of our body. This energises the body and contributes significantly to well-being. Most people occasionally suffer from impaired blood circulation, especially cold feet and hands are quite common. The back and neck warmer can help to improve your blood circulation.

Easy to use

through safe digital control

- 1. Place the back and neck warmer in position.**
- 2. Fasten it.**
- 3. Use the digital control to adjust the heat output and intensity of the massage function.**

Two heat settings

provide your body with pleasant warmth.

Press the left button to start the heat output.

Press the button again to increase the heat output.

Two massage intensities

You can choose between two massage intensities as per your individual preference. Press the right button to start the massage function. Press the button again to increase the intensity of the massage.

INSTRUCTION MANUAL

Please read these instructions carefully and save them for future reference.

Thank you for purchasing this wrap, the neck and back wrap that applies heat and massage to help relieve tension. Before using the product for the first time, please read all these instructions so you use the product correctly.

Warnings

Read the instructions carefully before using the device. Store them for future reference. It's absolutely vital for the user not to fall asleep while wearing the electric wrap. If the wrap is working for long time, the temperature selector should be set to the continuous position.

- Do not use the electric wrap in a folded position!
- Do not scrunch the wrap up!
- Do not stick pins into the electric wrap!
- Do not use the electric wrap on animals!

Electric Shock Hazard

Connect and use the device only if the mains supply complies with the information on device's rating plate. Check the electric wrap regularly to see if it shows signs of run out or damage. Do not connect the device to the mains if the power cable, temperature selector, plug or the wrap itself shows any visible signs of damage. If the device is broken, unplug it immediately. To prevent hazardous situations, any repairs or work on the device should be done by a qualified technician from the brand's official technical support service. Check the device's power cable regularly: make sure it never makes contact with hot objects or parts or close to other hot home appliance. Do not expose it to sharp edges or corners. Use the electric wrap only in dry rooms (Do not use it in the bathroom or similar place.) Do not use the electric wrap if it wet or damp! Ensure the device is completely dry before using it. The wrap creates low-intensity electrical and magnetic field that can have a negative effect on pacemaker functions. We therefore recommend contacting your doctor or the manufacturer of the electric wrap before using it.

Danger of Overheating

Do not use the electric wrap without the supervision of someone who knows how to use the device! Do not cover the temperature selector. To ensure the temperature selector works properly, it should never be placed on top of or under the electric wrap. Nor should it be covered by any other object. Never wrap the power cord around the temperature selector! Continuous use of the electric wrap at high temperatures may cause burns to the skins.

Danger of injury

The electric Wrap should not be used by young children. Do not allow persons with restricted physical sensory perception or mental abilities or with lack of experience and knowledge to operate the device unless they are supervised or have been instructed on the use of the appliance by somebody who is responsible for their safety. Keep out of reach of children. Supervise children to prevent them from playing with the device.

When storing the appliance, allow it to cool down before folding. Examine the appliance frequently for signs of wear or damage. If there are such signs, if the appliance has been misused or does not work, return it to the supplier before switching it on again. This appliance is not intended for medical use in hospitals. Children under the age of three are not to use this appliance due to their inability to react to overheating. The appliance is not to be used by young children over the age of three unless the controls have been pre-set by a parent or guardian, or unless the child has been adequately instructed on how to operate the controls safely. The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance. If the supply cord is damaged, it must be replaced by the manufacturer, its service agent or similarly qualified persons in order to avoid a hazard.

Never use the heat device if it is folded wrinkled, or pushed together.

Do not puncture the heating pad.

Only use the heating pad indoors.

Do not wash.

Double Insulation

This appliance is double insulated; therefore, no earth wire is required. Always check that the mains voltage corresponds to the voltage on the product rating plate.

The set includes

1x electric heating wrap
1x instruction manual

Technical Specifications

Voltage: 220-240V
Frequency: 50Hz
Power: 100W
Material: 100% polyester

Important

Do not use if wet!

Do not use with a helpless person, an infant or a person insensitive to heat!

Overlong application can result in skin burns!

Do not wash!

Do not wring!

Do not expose to the sun!

Clean the pad as per user manual!

Instruction for use

This electric heating wrap is specifically designed to fit around the neck, shoulders, and upper and lower back.

1. Remove the electric heating wrap from its protective plastic packaging and unfold it.
2. Place electric heating wrap on your shoulders and back. Make sure it fits around the neck.
3. Adjust the attaching systems at the neck and waist to ensure the appliance is in full contact with the body.
4. This appliance has two heat levels and four massage levels that can be used individually or in combination for customized therapy.
5. Plug the device into the mains. The red pilot light on the lower part of the temperature selector will come on when the device is plugged into the mains. To start the wrap, press the temperature button and/or the massage button to start their respective functions.
6. Press the temperature button to select the desired heat level. The green pilot light indicates low temperature and the red pilot light indicates high temperature.
7. Press the massage button to select the vibration intensity and vibration mode:
 - a. Low vibration-constant-green pilot light
 - b. Low vibration-intermittent-flashing green pilot light
 - c. High vibration-constant-red pilot light
 - d. High Vibration-intermittent-flashing red pilot

Note: the heat function and massage function can be operated individually or at the same time.

8. To turn off the wrap, press the temperature selector will start flashing if there is any kind of product error. If this happens, unplug the product from the mains and check the power cable the plug.

Burn Hazard

When wearing the electric wrap, it should not be too tight, but should not be too loose, but should be worn loose against the body. Never wear the device without outer cover.

Auto-Off function

This appliance will shut off automatically after two hours of continuous use of the temperature function and after 30 minutes of continuous use of the massage function.

Cleaning

Electric Shock hazard!

Before cleaning the wrap, unplug the power cable from the mains! IMPORTANT: Only the outer cover is washable. Do not wash the inner part of the wrap

The outer cover cannot be dry-clean.

Do not use bleach on the outer cover

Do not wring dry the cover

Do not get the inner case (connected to the power cable) damp or wet.

Unplug the wrap from the mains and remove the outer cover. Remove the cover by undoing the zip on the side.

The cover is machine washable in warm water. Only do air dry. When the outer cover is completely dry, inset the inner case inside to use the wrap again,

When not in use, store method/way/place:

- keep the wrap clean and dry.
- Store the device in a dry place when it is not in the use.
- Carefully fold the product to store it.
- Allow the appliance to cool before you fold
- Please do not use any chemical moth protection.
- Do not place any heavy objects on the device during storage, as this may damage the electric heating wire.

Disposal of used electronic devices

After its service life, please dispose of this device according to the directive on Waste Electrical and Electronic Equipment 2002/96/EC (WEEE).

The symbol on the product or its packaging informs you that this product cannot be treated as normal household waste, but must be returned at a collection point for recycling electronic and electrical devices.

Warranty and Repair terms

Please contact your dealer or the service centre in case of a claim under the warranty. If you have to return the unit, please enclose a copy of your receipt and state what the defect is.

The following warranty terms apply:

1. The warranty period for this products is 2 years from date of purchase. In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
2. Defect in material or workmanship will be removed free of charge within the warranty period.
3. Repairs under warranty do not extend the warranty period either for the unit or for the replacement parts.
4. The following is excluded under the warranty:
 - a. A damage which has arisen due to improper treatment, e.g. non-observance of the user instructions.
 - b. All damage which is due to repairs or tampering by the customer or unauthorised third parties.
 - c. Damage which arisen during transport from the manufacturer to the consumer or during transport to the service centre.
 - d. Accessories which are subjected to normal wear and tear.
5. Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim. If the supply cord is damaged, it must be replaced by the manufacturer, its service agent or similarly qualified persons in order to avoid the hazard.

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