**TRIGGERBOW®**

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**Instruction manual (Page 14)**

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**Manuale d'Istruzioni (<https://triggerbow.de/downloads/>)**

**Manual de instrucciones (<https://triggerbow.de/downloads/>)**

## Modell: TRIGGERBOW®

### Inhalt

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Model: TriggerBow®

## Content

1. Application notes .....	15
2. Self-responsibility.....	17
3. Instructions for use for operation.....	17
4. Designation and function .....	22
5. Cleaning and maintenance .....	22
6. Warning .....	23
7. Contact us .....	23
8. Warranty and Service.....	24

## Scope of delivery:

TriggerBow and instructions for use

**ENGLISH**

Read these operating instructions carefully and keep them for future reference. Please make the manual available to all users and follow all the instructions given in this document.

Dear customer,  
thank you very much for your purchase decision. With the TriggerBow, we deliver an outstanding quality product that you can enjoy for a long time.

**Areas of application (Wellness):**

- Massage treatment

**1. Application notes**

The TriggerBow is a massage treatment instrument. Use the product only for the back and neck massage. Please use the TriggerBow only for yourself. The product is intended for private use only. The method of treatment specified in the instructions for use must be observed.

The 22,000 mm stainless steel ball serves as a treatment piece. This is held by a magnet in the adapter device. This feature ensures the rolling function, but the ball can fall out of the holder in the event of lateral or excessive force.

The manufacturer is not liable for damage caused by improper use. If you are unsure whether the massager is suitable for you, consult your doctor.

### **WARNING**

*Do not use the TriggerBow ...*

- immediately after surgery,
- during pregnancy,
- If you have Injuries or infections in the body region to be massaged (e.g. open wounds, spondylarthrosis),
- for pain with unexplained cause,
- with limited perceptive faculty (e.g. pain-relieving medication, alcohol),
- for serious illnesses,
- for tumors in the treatment area,
- when taking anticoagulants,
- immediately after spinal surgery,
- for severe nerve diseases (e.g. sciatic pain) or acute inflammation.

**Consult your doctor if you are unsure.**

*Notes on safety:*

- Patients with cardiac pacemakers or implantable cardioverter defibrillators must not position the treatment piece in the region of the heart.
- Keep the packaging material away from children.
- Repair the device only if it is guaranteed to function properly afterward.
- Suitable from the age of 16.
- Application permitted from middle childhood (completion of 7 years of age) under adult supervision.
- Attention: Danger of suffocation due to small parts or packaging material.
- Use only when use is safe: Correct application and functional device.

## **2. Self-responsibility**

You decide on the pressure intensity during the treatment. This should be carried out in the comfort area. You can / should feel a gentle pain.

## **3. Instructions for use for operation**

Through selective pressure treatment, the tense back muscles are effectively treated. Of course, other muscles can also be treated with the treatment piece. It is best to free the upper body and apply it directly on the skin.

Muscle tension occurs faster than one would like. Unfortunately, these are not only an unpleasant, but above all a painful affair. Treat the pain points yourself:

By scanning the affected back muscles independently, you will find areas that are light to very painful and have increased tension. At these points (locally limited muscle hardening), for example, you carry out a pressure technique. They only use so much force that the perceived stimulus still has a pleasant character. The design of the TriggerBow is designed in such a way that punctual compression of the affected region is possible.

Maneuver the massager with both hands so that the treatment piece reaches the desired back section. The intensity of the treatment can be increased by leaning the device against a wall. Treat your own back muscles with pleasant massage techniques. These techniques have a relaxing effect. The device thus offers a punctual and powerful massage. An increased muscle tension is effectively reduced. Stop the massage if you find it painful or uncomfortable.

---

#### Pressure technique

---

Press on the locally limited muscle hardening and maintain this pressure for 30 - 60 seconds.

---

#### rolling motions


---

You massage the affected muscle area by rolling movements. It makes sense to combine the massage with circular movements. Do not treat the same body area longer than 10 minutes.

## Fields of application

Green: Treatment area

Red: Do not apply pressure / massage to fixed / bony areas.
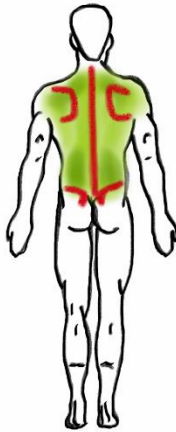
## Massage positions

The treatment is performed while sitting, standing or leaning against a wall. Breath calmly and try to relax.


You can treat yourself several times a day. After 2-3 days you should already notice the first improvements.

Reflect upon the results and adjust the treatment intensity accordingly.

The illustration below shows the **wrong application area!** Please do not press centrally on the spine, otherwise you may cause injuries.

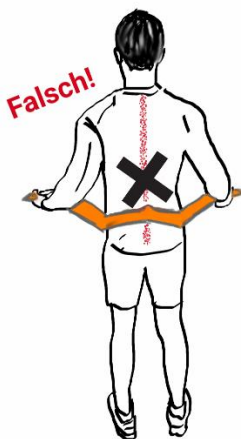There should be no disturbing or distracting factors during the treatment.

The TriggerBow cannot replace therapeutic treatment. If you have any questions, please contact your doctor.


*Unwanted side effects:*

Read this manual carefully. Pressure treatment should not be carried out on sensitive and bony regions of the body: In the middle of the spinal column (spinous processes), the iliac crest, sacrum or on the ribs.

Pay attention to an appropriate intensity, otherwise unwanted side effects can occur. With this device you treat exclusively the soft tissues (musculature) of the human body.

Improper use may result in bruising (hematomas), among other things. Make sure you use a sensible pressure during treatment so that you do not suffer any physical damage.

#### 4. Designation and function

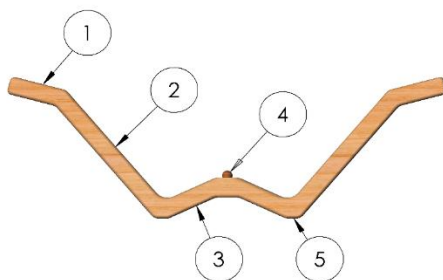- 
- 1, Grab handles:
  - 2, Long Arm
  - 3, Short Arm
  - 4, Treatment adapter: Ball + Magnet
  - 5, Anti-slip coating: Felt

#### 5. Cleaning and maintenance

Information on hygiene:

- Hot and wet can damage the wood surface
- Do not use chemical cleaners

Clean the treatment adapter with a damp cloth and dry it afterward. The device should not be stored in a humid environment. Store the device at room temperature to preserve its original shape. Any contact with damp surfaces can lead to unwanted changes in the shape of the product.

## 6. Warning

The product is made of wood and is therefore highly flammable. The lifetime of the TriggerBow depends on careful handling. Do not expose the device to strong temperature fluctuations, shocks or chemicals. Do not use abrasive cleaning agents. The cleaning process should be performed by an adult.

### **Disposal:**

Do not dispose of the device with your household waste at the end of its service life. Disposal can take place at appropriate collection points in your country.

## 7. Contact us

Please feel free to contact us if you have any questions or problems:

E-mail: [support@triggerbow.de](mailto:support@triggerbow.de)

TriggerBow

Streifing 14

94154 Neukirchen vorm Wald

[www.triggerbow.de](http://www.triggerbow.de)

### **Technology parameters:**

Name:	TriggerBow
Product type:	Massage device
EAN-Nr.:	4260594660011

## **8. Warranty and Service**

### **Quality certificate**

**We offer a 3-year warranty for material and manufacturing defects from the date of purchase.**

- Proof of purchase contract required
- Contact by phone, e-mail or Amazon
- E-mail with the subject "TB Garantie" to **support@triggerbow.de**

### **No warranty at all:**

- normal wear and tear  
(cracks and spalling, scratches)

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