
DermaWand[®] Pro

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Instructions for Use | Mode d'emploi |
Manuale di utilizzazione | Gebruiksaanwijzing

Gebrauchsanleitung	03
Instructions for Use	10
Mode d'emploi	17
Manuale di utilizzazione	24
Gebruiksaanwijzing	32

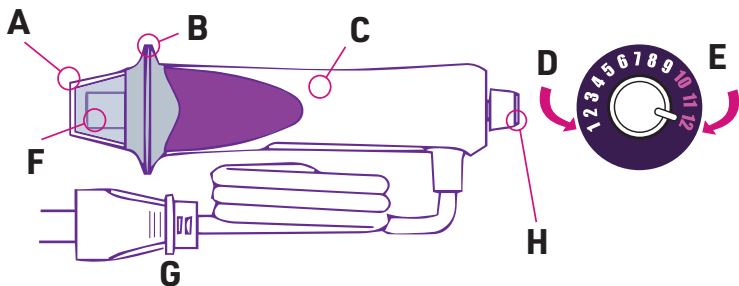DE

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IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

WARNING: READ ALL INSTRUCTIONS BEFORE USE

DANGER: To reduce the risk of electric shock:

- Always unplug the DermaWand® Pro from the mains outlet immediately after use.
- Never use the DermaWand® Pro while bathing.
- Never store or place the DermaWand® Pro in the bathtub or sink.
- Do not immerse or drop the DermaWand® Pro in water or any other liquid.
- Never attempt to retrieve a DermaWand® Pro that has fallen into water. Immediately disconnect the power plug.
- The DermaWand® Pro is a sealed appliance; do not attempt to open the casing. If it is damaged, call

Customer Services.

- The DermaWand® Pro should never be left unattended when plugged into a mains outlet.
- Extreme caution is advised when the DermaWand® Pro is used in the proximity of children or people with special needs.
- Only ever use the DermaWand® Pro for the purposes described in this manual. Do not use any attachments other than those recommended by the manufacturer.
- Never operate the DermaWand® Pro if the power cord or plug is damaged, if it does not work properly, if it has been dropped or damaged, or if it has fallen into water. Make sure to send your DermaWand® Pro to a service centre for inspection and repair by calling Customer Services.
- Keep the cable away from heated surfaces.
- Never use the DermaWand® Pro when you are overly tired or feeling unwell.
- Never drop or insert any object into an opening of the DermaWand® Pro.
- Do not use DermaWand® Pro outdoors, in areas where aerosol products (sprays) are used, or in which there are elevated oxygen levels.
- Do not use the DermaWand® Pro near flammable liquids or materials.
- Only ever connect the DermaWand® Pro to a properly earthed power outlet. Do not use the DermaWand® Pro if you are pregnant or planning to become pregnant, if you have a pacemaker, or if you have varicose veins or any other skin condition.
- Do not use the DermaWand® Pro on mucous membranes (eyes, ears, nose, mouth).
- Never use the DermaWand® Pro internally.
- Dry your hands completely before using the DermaWand® Pro.
- Never hold the DermaWand® Pro by the toning tip. Always keep your hand behind the connection.
- Do not use the DermaWand® Pro if the toning tip is broken.
- This appliance may be only used by children from 8 years of age or by persons with reduced physical, sensory or mental capabilities or lack of experience and/or knowledge if they are supervised or have been

instructed in the safe use of the appliance, and can understand the hazards that may arise from its use. Children should not play with the appliance. Cleaning and user maintenance should not be performed by children without supervision.

RETAIN THESE INSTRUCTIONS.

The DermaWand® Pro is a hand-held device that emits an electrostatic current. This appliance is not a prescribed treatment, and is intended for household use only.

During operation, the DermaWand® Pro continuously emits electrostatic impulses to the skin for as long as the toning tip touches the skin surface. The knob positioned at the rear of the device controls the intensity of the impulse delivered to the skin surface.

Certain exercises can be performed with the device, with the toning tip being brought into contact with the skin, passed over the skin surface in a particular area, and then held for a short period of time before being raised from the skin surface.

The potential result of regular, routine use of the DermaWand® Pro is a reduction in the number and size of fine lines and wrinkles, which in turn can improve the overall appearance of the skin.

- **FEELING:** When the DermaWand® Pro stimulates and massages your skin, you will feel a tingling sensation that increases as you intensify the toning strength. Follow the instructions contained in the manual carefully, in order to ensure that you select the intensity that suits you best.
- **NOISE:** During your treatments, you will hear a deep humming-like sound.
- **ODOUR:** The DermaWand® Pro provides instant stimulation and heat massage. This process releases enriched oxygen, which creates a fresh natural odour.
- **OPTICS:** When using your DermaWand® Pro, its tip will glow. There may be a very slight redness of the skin; this is normal. It is very **IMPORTANT** to apply a good topical cream **BEFORE** using your DermaWand® Pro. If more noticeable reddening of the skin occurs, discontinue use of the DermaWand® Pro for a few days until the redness disappears. Start again with a lower setting and ensure that you apply a good moisturiser both before and after treatment. Do not use products containing fruit acids, glycolic acid, or other abrasive products.
(Note: Discontinue use if excessive redness or swelling occurs.)

The results vary from person to person, depending on age and the length and frequency of treatments. You should

be able to see and feel results within a week!

Before using your DermaWand® Pro

- Before plugging in your DermaWand Pro, REMOVE the protective cap from the toning tip, and inspect it for cracks or damage.
- Turn the toning intensity down to the lowest setting level. We recommend that you start with the lowest setting for the first few days, until you get used to the feeling. When using the device for the first time, first try it out on the back of the hand to accustom yourself to the sensation.
- Plug the DermaWand® Pro into a suitable power outlet and allow the device to warm up for one minute before use. You will notice an orange glow in the toning tip casing.
- Be sure to hold the appliance with your fingers on the grips, and make sure that your fingertips never go in front of the connector.

After treatment with your DermaWand® Pro

REMOVE THE PLUG FROM THE MAINS BEFORE replacing the cap for storage.

Getting Started

It is recommended that you perform treatments with your DermaWand® Pro in two phases.

PHASE 1 is the 3-minute swipe. This phase generally covers the entire face. It is recommended that this phase be completed at least twice a day.

PHASE 2 is the facial exercise routines. These exercise routines are tailored to specific areas of your face.

PHASE 1: 3-minute swiping treatment
Follow these steps twice a day (morning and evening).

STEP 1: When using the DermaWand® Pro for the first time, turn the knob to level 1. Be sure to remove the safety cap first. First try it on the back of your hand in order to trial the tingling sensation at the selected intensity level. If you have never used your DermaWand® Pro before, or if your DermaWand® Pro has not been used for several days, it may take up to 1 minute for the toning tip to begin to glow.

STEP 2: Apply the toning tip directly to your facial skin (remember to remove the cap). For the first few days always apply with the lowest setting, until you get used to the feeling. If you no longer feel the intensity of the stimulus (a slight tingling sensation), it means that your skin is ready for the next stage. Increase the intensity on the selector knob a little further, until you reach the limit of your comfort zone. The DermaWand® Pro is effective even at the lowest setting, and your comfort level depends on the sensitivity of your skin.

STEP 3: Start by gently swiping the DermaWand® Pro over the entire face and neck for three minutes in a straight or circular motion. Be sure to keep the DermaWand® Pro in constant motion across your face. Focus on areas of very deep wrinkles for 30-40 seconds each. Lines and wrinkles take longer to respond.

STEP 4: When you are finished with the general facial swipes, start with some of the facial exercises shown in this guide. The more exercises you perform each day, the faster you will see results. If reddening of the skin occurs, discontinue use of the DermaWand® Pro for a few days until the redness disappears. Start again with a lower setting and ensure that you apply a moisturiser before using the appliance. Also apply a moisturizer or facial cream after using the DermaWand® Pro.

CAUTION: If excessive redness or swelling occurs, discontinue use immediately.

PHASE 2: Facial exercise routines

For maximum results, perform the following exercises twice a day.

NECK: Apply a moisturiser, then beginning on the neck area, slowly slide the device up to the chin. Then return to the neck area and repeat 6-8 times, performing a full ,V-shaped' movement.

CHIN AND UNDERCHIN: Apply a moisturiser to the chin area. Now start in the middle directly under the chin, and move the device outward, along the chin arch and up to the earlobe in small, tight circular motions, in a circle approximately the size of a one-pence coin. Do this for 30-45 seconds on each side. In addition, you can treat this area with a second exercise, starting in the middle of the chin and sliding to the base of the ear, holding for 3 seconds, then repeating 2-3 times on each side. Always apply a good cream or lotion before using your DermaWand® Pro. Generally use a moisturiser for the face and forehead, and an eye cream for the

area around the eyes, including any crow's feet.

WRINKLES AROUND THE MOUTH:


Apply a moisturiser and then begin on one side at the lower centre of the chin. Move the device in an upwards direction with a zig-zag movement, just above the corner of the mouth. Repeat this exercise for 30 seconds, then perform it on the other side for 30 seconds as well.

NASOLABIAL FOLDS: Apply a moisturiser and then perform zig-zag motions with the device, starting from the bottom of the fold, and following the profile upwards. Perform this exercise on one side of the face for 30-40 seconds, then repeat on the other side of the face.

CHEEKS: Apply a moisturiser and move the device from the joining point of the nose in a circular motion across the

cheek to the anterior edge of the ear. Perform this exercise for 30 seconds, then repeat on the other side of the face.

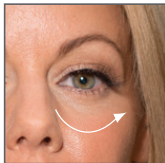LINES UNDER THE EYES: Apply a good eye cream before this exercise. Move the device from the inner corner of the eye to the outer corner of the eye in a single sliding motion while pressing the device to the bone area. Perform this exercise for 30-40 seconds, then repeat around the other eye.


DO NOT SET THE TONE CONTROL KNOB ANY HIGHER THAN LEVEL 3 FOR THIS EXERCISE.

CROW'S FEET: Apply good eye cream before this exercise. Perform circular motions across this area. Do this for 30-45 seconds on each side.

EYEBROW AREA: Apply a moisturiser, then slowly move upwards across forehead from the inner edge of the eyebrow to the hairline and hold for three seconds. Repeat this movement from the middle of the eyebrow to the hairline, and then from the outer edge of the eyebrow to the hairline. Repeat this three-step exercise 5-6 times on each eye. You may notice a slight redness the first few times you perform this exercise.

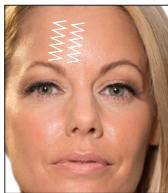FOREHEAD: Apply a moisturiser and, with moderate pressure, apply the DermaWand® Pro, pulling the skin of the forehead in a slight zig-zag motion toward the hairline. Once you have reached the hairline, hold the skin taut with the appliance for 15 seconds. Perform this exercise on each side of the face three times, starting from the

inner edge of the eyebrow, then from the middle of the eyebrow, and finally from the outer edge of the eyebrow.

BETWEEN THE EYEBROWS: Apply a moisturiser and move upwards in a zig-zagging motion. Repeat this exercise for 45 seconds.

Broken toning tip

The toning tip is not covered by the DermaWand® Pro's product warranty. Always place the safety cap on the toning tip when not in use. Be careful not to drop the DermaWand® Pro during use, as the toning tip may break.

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EN. Do not dispose this product at the end of its service life through normal domestic waste. It has to be disposed through a collection agency that recycles electric and electronic devices. This is pointed out by the symbol on the product, in the instruction manual, and on the packaging. Please inquire about local collection agencies that may be operated by your distributor or your municipal administration. By recycling, making use of the materials, or other forms of recycling old units you are making an important contribution to protecting our environment.

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Hergestellt in China | Made in China | Fabriqué en Chine | Fabricato in Cina | Made in China

**Importeur:**

CH: MediaShop AG | Industriering 3 | 9491 Ruggell | Liechtenstein
EU: MediaShop GmbH | Schneiderstraße 1, Top 1 | 2620 Neunkirchen | Austria

DE | AT | CH: 0800 376 36 06 – Kostenlose Servicehotline

ROW: +423 388 18 00

office@mediashop-group.com | www.mediashop.tv