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Shiatsu foot massage device

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Shiatsu foot massage device

Features

- Combination of air pressure, Shiatsu, foot reflexology, kneading massage with integrated heat function
- Air massage intensities 6 levels adjustable
- Massages the entire foot up to the toes
- two heat levels adjustable
- automatic power-off
- manual and automatic program selection
- safe power supply via DC12V adapter
- Portable and practical

Before first use:

Unpacking

Remove the device and the operating instructions from the box.

Before using the device for the first time, remove all packaging materials and protective films from the device, especially the protective film on the underside of the device.

NOTE

Keep the original packaging during the warranty or guarantee period of the device in order to be able to send the device properly packed in the event of a warranty or guarantee claim. Transport damage will invalidate the warranty or guarantee claim.


Check the scope of delivery

Check the contents for completeness and for visible damage. Report an incomplete or damaged delivery to your supplier immediately.

1. DC socket
2. Power adapter
3. DC connector
4. Display screen & control panel
5. Airbag position
6. Removable cover
7. Venting the security key

Control Panel Instruction

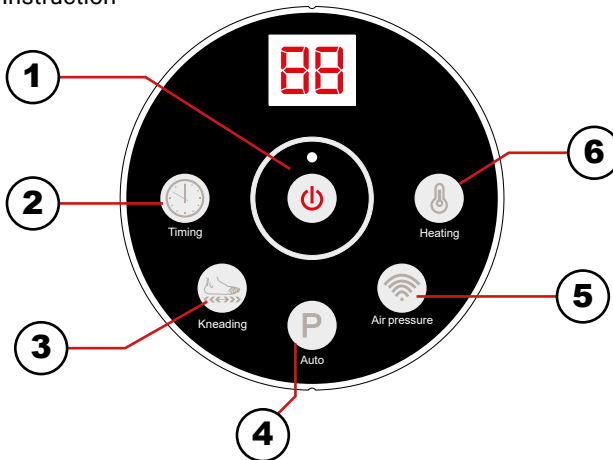1. Power switch: On/Off
2. Time setting button: Selection of the time setting
3. Kneading button: Switch on or off for kneading massage
4. Auto key: Switch on/off | automatic mode
5. Air pressure button: Set / switch air pressure level on and off
6. Heat controller: Low or high temperature can be selected manually or automatically

Commissioning and use:

Connect the supplied power cable / 12 volt adapter to the connection of the foot massager and the power plug of the supplied power cable / 12 volt adapter to the power outlet. A signal tone sounds and an LED lights up on the control unit. The device is now ready for use.

When switching on the power switch, key „1“, you are initially in manual setting mode.

With key „2“ you can determine the duration of the selected program. (10min /20min /30min for selection).

With key „5“ you can regulate the air pressure intensity. You can choose between 6 levels.

With key „6“ you have the possibility to switch on the heat function if desired. This is possible in two steps (low/high). Of course, this function can also remain switched off.

With key „4“, labelled „P“, you can switch between 2 automatic programs.

P1: Combined air massage mode with kneading massage mode/air pressure massage/kneading massage, the three above mentioned function types in the cycle

P2: Combined air pressure massage/kneading massage/air massage mode with kneading massage mode, the above three types of functions in the cycle

Press key „4“ again to switch off the automatic cycle.

When the massage time up to the setting time is reached, the buzzer emits the tone „di“ three times and the device stops working.

Product maintenance and safety instructions:

In order to keep the massager in the best condition, please make sure it is used correctly and observe the following points.

NOTE

Never massage feet with injuries or illnesses of any kind. Use the foot massager only when seated.

Please do not use the massage continuously for a long period of time at the beginning. You could increase the application time after you have become accustomed to it. If abnormal conditions such as overheating or odour occur during use, please stop using the foot massager immediately and switch it off. After use, please turn off the system by first pressing the „Button „1“ on the control unit and then unplugging the power cord from the wall outlet. Do not pull out the power cord by force.

Please unplug the power cord from the power outlet before cleaning the product. Please use a soft cotton cloth and a neutral detergent to clean the massager. If you are under medical treatment or have any of the following, please contact your doctor and follow your doctor's instructions.

- Cancer
- Chronic or serious illness
- Severe heart disease
- Pacemaker
- Other metal utensils, such as prostheses
- Mental illness
- Pregnant women or children

Please do not pull the power cord of the massager by force.

Please do not use the foot massager in water or damp places. If the foot massager accidentally falls into the water, please be sure to turn off the power and unplug the power cord before removing the massager from the water. Please do not place the product under a pillow or blanket when using it. If you cover the massager, the motor temperature is too high and may cause malfunction or damage. Please do not use the massager near ovens or other heat sources. Please do not unscrew the foot massager yourself if it is defective. Such behaviour will result in the loss of your warranty. For your personal safety, please do not unscrew the zipper. We do not take responsibility for the consequences of such behaviour.

NOTE

Basic Safety Instructions

Check the device for visible external damage before use. Do not operate a damaged or dropped device.

Persons who are not able to operate the device safely due to their physical, mental or motor skills may only use the device under the supervision or instruction of a responsible person

Allow children to use the device only under supervision.

Have repairs to the appliance carried out only by authorised specialist companies or customer service. Improper repairs can cause considerable danger to the user. In addition, the warranty will become void.

Protect the device from moisture and the penetration of liquids. Avoid direct sunlight.

NOTE

Requirements for the installation site

For safe and fault-free operation of the unit, the installation site must meet the following requirements:

The unit must be set up on a firm, flat and horizontal surface. Do not set up the device in a hot, wet or very humid environment or near inflammable material.

Conversion or alteration of the product will impair product safety.

Caution! Risk of injury! Never open the product unauthorized. Never carry out repairs yourself! Do not operate under cover.

Handle the product carefully. It can be damaged by shocks, blows or falling from a low height. Do not operate the product once there is visible damage. Keep the product away from moisture and extreme heat. Never immerse the product in water or other liquids

Warranty/warranty

The warranty period is 24 months from the date of purchase. This includes a 6-month guarantee from the date of purchase.

The territorial scope of the warranty protection is Europe-wide.

The proof of purchase is valid as proof of this. The product has been thoroughly checked for quality and functionality before dispatch.

The guarantee is excluded in the following cases:

- in the event of modifications to the system without our approval
- in the case of use other than that stated above
- in the event of damage caused by falling

When making a claim under the guarantee/warranty, please return the complete product in transport-safe packaging to the point of sale from which the product was purchased. Describe the detected defect in the product and enclose the original proof of purchase with the shipment.

The seller reserves the right to remedy the defect before refunding the purchase price.

Necessary transport costs will be reimbursed to the customer if a justified claim under the statutory guarantee/warranty is established.

The use of the legal rights for defects is free of charge. These are not restricted by the guarantee.

Mobiset GmbH hereby declares that this product complies with the essential requirements of the applicable European directives.

Intended use

This device is only intended for use as a foot massager.

Any other or additional use is considered improper.

Simple malfunction analysis and how to deal with it

Phenomenon	Possible Cause	Possible help
Not working?	1. power supply is not switched on	1. Connect power supply
	2. mains plug has bad contact	2. check plug
	3. the "Power" key was not pressed	3. press the „Power“ key
Stopped working	1. switched off	1. check if the power supply is secured
	2. the massage time is over 15mins	2. press the „Power“ button again

Specification

Product name: Seecode Shiatsu foot massage device

voltage: 12V

Power: 40W

ATTENTION:

1. This device is not intended to help people with physical disabilities, sensory or nerve

defects, or lack of experience and common sense to use this device unless they have been supervised and advised on how to use this device by the persons responsible for their safety.

- Children must be supervised to ensure that they do not take any instruments to play with.

Care instructions for the removable footmuff by means of a zip

Wash at a temperature of up to 60° C in a gentle or normal wash cycle, e.g. a delicate wash or easy-care programme.

The footmuff must not be treated with bleach, it should only be washed with colour detergent and mild detergent.

Dry exclusively in the easy-care programme of the dryer.

The footmuff must not be ironed.

Dry cleaning is not suitable for this product.

Further information

The foot reflexology massager is made of a shatterproof ABS plastic on the outside and the inner bag is made of a removable and washable fabric cover.

Care instructions foot reflexology massager

Please do not use any bleaching, caustic or abrasive cleaning agents. Hot water or a simple washing-up liquid is sufficient to remove any dirt. It is recommended to wash the footmuff inside out.

Removing the footmuff

The footmuff is attached to the fixed fabric frame by a zip. The footmuff is also fixed inside the unit by a Velcro fastener. After opening the zip, please remove the footmuff. When reinserting the footmuff, please fix the Velcro fastener inside and close the zip fastener again.

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