Digitaler Mundgeruch Tester

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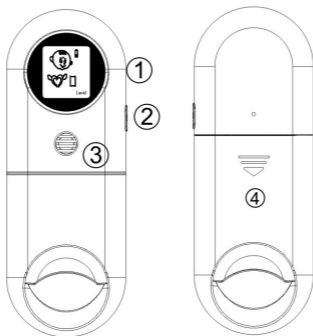
Digital Pocket Bad Breath Tester Breathalyzer

(1) LCD display
(2) Power button

(3) Breath Inhaler
(4) Battery compartment

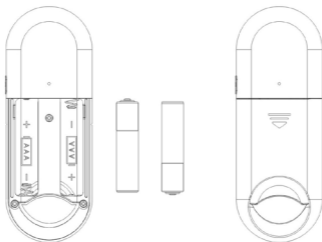
Specifications

- Sensor type: Semiconductor sensor
- Warm up time: 9 S
- Blowing time: 5S
- Working Environment: Relative Humidity:
≤95% no dews
- Indication manner:
Digital and Cartoon LCD indication with backlight
- Power: 2x AAA (batteries not included)
- Dimension: 115*45*15mm; net weight: about 28g

Cautions

The following situations could affect your result. Use the device after the displayed time below.

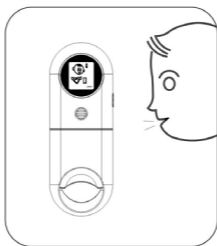
Brushing teeth, mouthwash	Within 30 minutes
Consuming candy or Gum	Within 30 minutes
Smoking	Within 30 minutes
Food	Within 30 minutes
Drinking or the morning after	Within 24 hours
Coffee, tea, or beverage other than water	Within 30 minutes

1. Insert two "AAA"(UM-4)alkaline batteries with correct polarity as indicator.
2. Replace the battery compartment lid.

Warm Up

1. Press and hold the power button for about 2 Seconds. the LCD will turn on with one beep sound.

When the screen shows “c” symbol now user blow to gas entrance about 5 seconds.

2. If the first time warm up is failed the tester will repeat warm up till the sensor is cleaned.

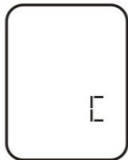
Note

It always has a large number of gas and impurity adsorbed on the surface of odor sensor, when first time taking out the tester from the packaging, or after long time no using for the tester, it need more time for cleaning and warm up the sensor, you should repress the power button and retry 3-5 times if possible.

If you blow and get 3-4 results, should be repress the power button and retry 3 time till the sensor is cleaned.

Test

1. The “C ”symbol is displayed with one beep sound when warm-up is succeed, and the 9 to c countdown is started that indicator the tester is in warm-up mode.

2. Near and breath softly into the breath inhaler for several seconds, suggest to blow about 5 seconds.
3. After one beep sound, look into the test result from the display.

Very good

Normal

Not so good

Bad

Very bad

Operation

- 1.1. Open its back battery compartment, put the batteries according to its polarity then close down it.
- 1.2. Press On/Off button and hold for 2 seconds. the tester enter into warm-up status, the screen will show 9, 8. C count down.
- 1.3. When warm-up finished, the screen will show "C Twinkling, now user should blow to gas entrance about

5 Seconds, we recommend the blowing distance is between 5-15mm.

1.4. The tester will give out a "di" sound, then it will shown testing result will keep on 3 seconds

then power off automatically.

1.5. If new testing needed, please turn on again

Notification

1. Avoid any strong shock and fall

2. This tester is sensitive to some polluted factors. Keep the unit away from all liquids and store in a safe environment, if noise gas exist, the normal testing would be disturbed.

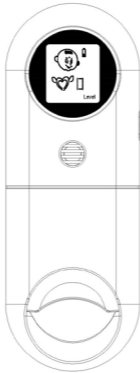
3. Wait at least 10minutes to test after cleaning your teeth or using mouth wash to keep its accuracy

4. If the tester was not used for a long time, try to test it several times before normal testing to re-activate it

Fixing the hand strap

1. Insert the smaller loop into the hole in the housing

2. Thread the longer loop through the eyelet and pull through

Backlight

The tester have no back light in default.

- Portable design
- Healthy No-touch design
- Battery saved design, low battery indication
- Advanced Nanometer semiconductor gas sensor
- Digital LCD and Cartoon indication with backup light
- Quick response and resume
- Self-checking to show different grad

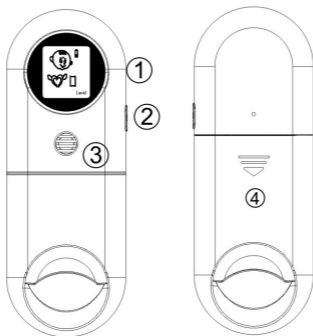
Features:

Everybody worries that bad breath would make people feel uncomfortable, especially girlfriends or boyfriends when face to face together, they are sensitive of your breath, but how to avoid these Embarrassing scene???

Try this Odor Breath Tester!!!! It will help you to know your breath condition is good or bad!!!!

Odor Breath tester is specially designed to check people's breath quality. It adopts semiconductor gas sensor, which with excellent sensitivity and responsibility, simple and convenient to use. It will show the final testing result in digital and cartoon mode, everyone can check themselves by this tester.

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Mobiset GmbH
Mobiset GmbH - Rösrather Straße 333 - 51107 Köln