HERO **FX**

HIGH FREQUENCY MASSAGE DEVICE

USER MANUAL

ACTIVE  HERO

Summary

Get started with HERO FX	3
Packaging content.....	5
Description of the HERO FX.....	6
Getting started	7
Instructions for use.....	9
Different massage heads.....	10
Different vibration speeds.....	10
Areas of use of the HERO FX.....	11
Advice on care, safety and tips for use.....	12
Technical features	14
Disposal.....	14

Massage guns are now an integral form of sports reinforcement exercises and are used by sports professionals, weight lifters and endurance athletes.

Get started with HERO FX

Physiotherapists, chiropractors and sports coaches are turning to vibration therapy to help athletes improve their performance and recover from their injuries. **HERO FX** is your new fitness ally.

Use the **HERO FX** before your workouts to increase your performance or after your workouts to speed up recovery. The vibration technology of the **HERO FX** penetrates deep into your muscles and fascias.

The high frequency pulsations of **HERO FX** can help relax contracted muscles, increase blood circulation in the area and work on your fascias.

The **HERO FX** is easy to carry and its lithium-ion battery offers a range of up to 3 hours. Different massage speeds and pressure intensities allow you to select a power adapted to each of your needs.

Give your training a new dimension and discover the power of HERO FX!

Go through the exercises described in this booklet to find out how to get started with your **HERO FX**.

ACTIVE  HERO

HERO FX

HIGH FREQUENCY MASSAGE DEVICE

Packaging content:

- HERO FX – High frequency Massage Device
- 2 x batteries
- charger
- transport case
- 5 x massage heads
- instructions for use

Please read this user manual carefully, keep it for future use, make it available to other users and follow its instructions. Follow the steps in the **Getting started** section on page 7 to mount and load your HERO FX.

Description of the HERO FX

LCD display screen

Displays current running time, RPM, applied force, and battery display for 4 levels

- Button

Reduce vibration speed

ON / OFF Switch

Press the ON/OFF switch to light up the HERO FX

+ Button

Increasing vibration speed

Charging point to connect to the battery power supply unit

Massage head (1 of 5)

LCD display screen and buttons

Battery

Handle

Getting started

Fully charge the battery

- Fully charge the battery before using the **HERO FX** for the first time (up to 6 hours)
- To do this, connect the power supply to the charging socket of the battery
- Plug the adapter into a power outlet
- When the battery is fully charged (the power supply LED turns from red to green), unplug the power supply and firmly insert the battery into the **HERO FX** handle

Attaching a massage head

- Once your **HERO FX** is fully charged, you can attach a massage head
- Choose the preferred massage head
- Place the massage head so that the small nodule on the massage head is facing the notch on the massage head on main part of the **HERO FX**
- Push the massage head firmly in to fix it in place

Also see the different **massage heads and vibration speeds** on page 10.

Instructions for use

Read the instructions before using the **HERO FX**. Keep this manual for future reference.

Charging

- Completely charge the **HERO FX** before its first use (up to 6 hours)
- After the first charge, the **HERO FX** will require no more than around 1 hour and 30 minutes charging time
- The battery is fully charged when the LED on the power supply goes from red to green
- Maximum usage time extends up to around 3 hours (depending on the speed used)

Using HERO FX

- Press the ON/OFF button to turn on the **HERO FX** at speed level 1.
Select from one of the 15 speed vibration levels:
 1. Press the + button to increase the vibration speed
 2. Press the - button to reduce the vibration speed
- Move the massage head along a muscle group at the desired speed.
Use different levels of pressure according to your needs

Removing and replacing a massage head

- Press the ON/OFF button to turn off the **HERO FX**
- Grasp the attached massage head and pull firmly to remove it
- Insert another massage head in the device
- Be sure to place the small nodule of the massage head facing the notch on the massage head of the device
- Push the massage head firmly inside to secure it

TIPS FOR USING THE HERO FX

- Move the **HERO FX** over the desired muscle groups for between 15 seconds and two minutes
- Adapt the pressure applied according to your needs, but never use excessive force
- Never use the massage head for more than 3 minutes on the same area!
- Prolonged use could cause tissue damage

Different massage heads

Five interchangeable massage heads are included with your **HERO FX**. These massage heads are designed to treat different areas of your body and meet specific needs for therapeutic massage.

Round ball, small

For recovery and reduction of body aches after training.
Ideal for small muscle groups.

Round ball, big

For recovery and reduction of body aches after training.
Ideal for large muscle groups.

Cone

Specifically designed for problem areas and the treatment of « trigger points »

Flat

For muscles and areas that need special attention.
Simulates a deep massage.

Spinal


Ideal for back muscles situated along the spinal column.

Different vibration speeds


The 15 different massage speeds make it possible to select the right intensity for a specific muscle group or trigger point. Each speed offers a different vibration frequency (Hz) and a different number of revolutions (rpm).

The number of revolutions can vary between 1.600 and 3.000 rpm.


Areas of use for the HERO FX

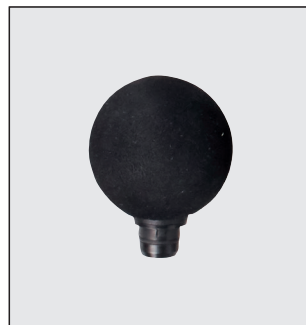**For big
muscle groups**

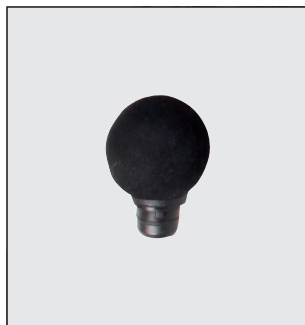
**For small
muscle groups**

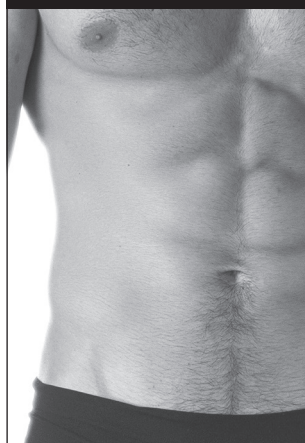
**For all parts
of the body**

**For punctual
treatment**

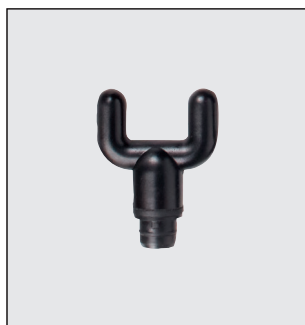For back muscles

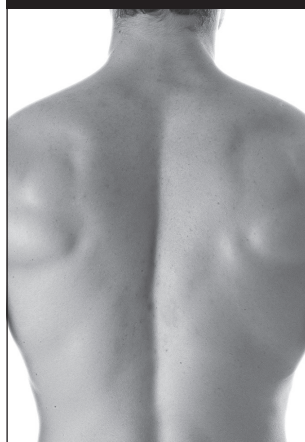Advice on care, safety and tips for use

Maintenance, cleaning and storage

- Use a soft cloth to wipe off dirt or moisture on the device
- When you want to store or transport it, remove the battery and store it in the included case

Safety notes

- Please read the following warnings carefully to minimize the risk of burns, fire, electric shock or injury
- **Seek medical advice before using this product. Do not use it if you are injured**
Exclusively for adults. Keep out of the reach of children and pets
- **The use of this device is not a substitute for a medical consultation or treatment.**
In case of any pain or illness, consult your doctor immediately!
- Use the **HERO FX** exclusively according to the instructions provided in this manual
- Keep away from liquids and heat sources
- Never let the **HERO FX** run unattended
- Unplug the appliance from the mains after charging or before use
- Do not remove any screws or attempt to disassemble the device
- Do not exceed one hour of consecutive operation. After one hour of operation, leave the unit off for 30 minutes before using it again

For the sake of your health, it is not recommended that you use the massage gun in the following cases:

- Presence of electrical implants (such as a pacemaker, for example)
- Presence of metal implants
- Whilst wearing an insulin pump
- In case of fever
- In case of known or acute heart rhythm disorders and other disorders of cardiac excitability and conduction
- In case of seizures (eg epilepsy)
- Pregnant women
- People suffering from cancer
- Following operations where increased muscle contractions may interfere with recovery
- Never use the device near the heart. Stimulation electrodes should never be placed in front of the ribcage (the ribs and sternum of which form the visible part), and especially not on the two large pectoral muscles. Such use may increase the risk of ventricular fibrillation and lead to cardiac arrest
- On the skull, in the mouth, throat or larynx
- On the neck / carotid
- On the genitals
- On skin with severe or chronic lesions or inflammations, eg. in case of inflammation, redness, rash, burns, bruises, swelling, painful or not, and in case of open wound or healing, on postoperative scars during cicatrisation
- In environments with high humidity, eg. in a bathroom, in the bath or in the shower
- Do not use after consuming alcohol

Important precautions for use

- Only use **HERO FX** on soft tissues
- Do not use **HERO FX** on your head or other bone areas
- Hematomas may appear, regardless of the chosen rate of vibrations
- Control the areas you treat frequently and adapt your use at the first signs of pain or discomfort
- Use only on clothes or clean dry skin
- We recommend that you use the **HERO FX** over your clothes (eg. tight fitting sportswear, leggings, etc.)
- Keep your hair, fingers, and other body parts away from the base of the massage head to prevent bruising

WARNING!

Use this device exclusively for massage:

- On people
- For the purpose for which it was developed and in the manner indicated in this instruction manual
- Incorrect use can present risks
- External use
- With the original accessories supplied and re-orderable, otherwise the right to the guarantee will not apply

Technical data

Article number:	21389 / HERO01
Model:	WO2
Dimensions:	Ca. 20 cm x 17.5 cm x 10.5 cm
Battery:	Lithium-Ion
Battery lifetime:	Up to 3 hours (depending on the speed used)
Battery charging time:	Around 90 Minuten
Power supply:	220 - 240 V ~ 55 Hz
Power:	1500 mAh 12 V
Speed:	15 speed levels up to 3000 Rpm
Noise level:	Approximately 55 DB at maximum Speed
Weight:	Approximately 900 g (battery incl.)
protection class:	II
Made in PRC	

Disposal

Packaging materials are recyclable. Dispose of the packaging in an environmentally friendly way and take it to a recycling centre.

Dispose of the device in an environmentally friendly manner. It should not be thrown away with household rubbish. Dispose of this unit in a recycling centre for used electrical and electronic equipment. For more information, consult the local guidelines of your borough.

Broszio & Co.
GmbH & Co. KG
Ottensener Str. 1-5
D-22525 Hamburg
Tel.: 040 - 5472920
E-Mail: info@broszio.com
www.broszio.com

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